



Horizontal Electrophoresis

DS-SHET-B Horizontal Electrophoresis Tank

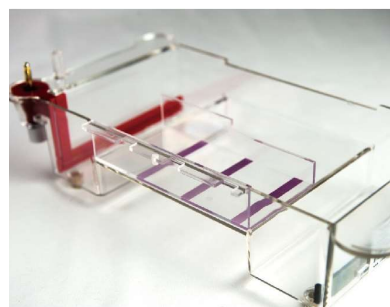

With the drift prevention function on the tray, effectively prevent the band bending.

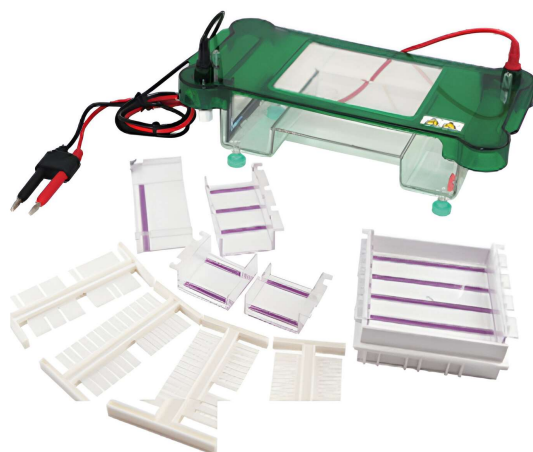

Product Features

Abundant buffer not only guarantees cooling effects, but also keeps the pH value stable during the whole experiment process.

Removable electrodes enable maintained and replacement easy and quick.

With background color design on the tray, convenient to adding samples.

It's more convenient to open the lid with the easy-open buttons.

It prevents buffer volatilizing and electricity leakage with transparent lid.

Auto-switch-off when the lid is moved.

DS-SHET-B Technical Specifications

MODEL	DS-SHET-B
Gel size (WxL)	120x120 (mm)
	120x60(mm)
	60x120 (mm)
	60x60 (mm)
No. of samples	1~100
Total buffer volume	~550 (ml)
Dimension (LxW xH)	300x160x75 (mm)
Weight	1.1 (kg)

DS-SHET-B Packing List

Parts name	Quantity
Lid with cables	1
Body tank with electrodes	1
120x120mm Gel tray	1
60x120mm Gel tray	1
60x60mm Geltray	1
Gel box	2
1.5mm thickness 1/2/3 wells comb	1
1.5mm thickness 1/2/3 wells comb	1
1.0mm thickness 8/18 wells comb	2
1.0mm thickness 1/25 wells comb	2
level adjusting button	2
level adjusting button	4