

# Horizontal Electrophoresis

## DS-SHET-F Horizontal Electrophoresis Tank



### Product Features

Abundant buffer not only guarantees cooling effects, but also keeps the pH value stable during the whole experiment process.

Removable electrodes enable maintenance and replacement easy and quick.


Adding sample more quickly and suitably for the multi-channel pipettes.

With a special gel casting, preventing the gel leakage, operation is simple and convenient.

It's more convenient to open the lid with the easy-open buttons.

It prevents buffer volatilizing and electricity leakage with transparent lid.

Auto-switch-off when the lid is moved.

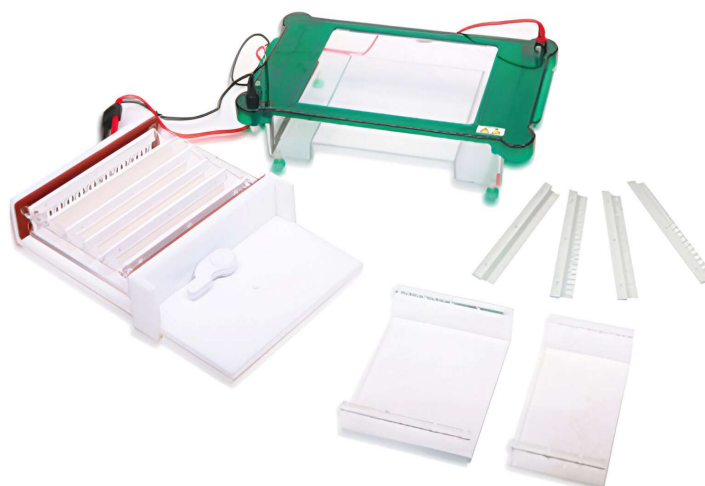

### Product Overview

#### DS-SHET-F Technical Specifications

MODEL	DS-SHET-F
Gel size (WxL)	200x200 (mm)
	200x150 (mm)
	200x100 (mm)
No. of samples	17~104
Total buffer volume	~800 (ml)
Dimension (LxW xH)	397x230x93 (mm)
Weight	3.3 (kg)

#### DS-SHET-F Packing List

Parts name	Quantity
Lid with cables	1
Body tank with electrodes	1
200x200mm Gel tray	1
200x150mm Gel tray	1
200x100mm Geltray	1
Gel casting stand	1
1.0mm thickness 17-wells comb	2
1.5mm thickness 17-wells comb	2
1.0mm thickness 22-wells comb	2
1.0mm thickness 22 wells comb	2
1.0mm thickness 36 wells comb	2
1.0mm thickness 30 wells comb	2
1.5mm thickness 44-wells comb	2
1.0mm thickness 44-wells comb	2
Level adjusting button	4