

Horizontal Electrophoresis Tank



India Office

16 S/F B/S, T-Block Extn.
Jain Colony Part-2 Uttam Nagar,
New Delhi-110059
info@darsunscientific.in
(+91) 9999136670, 7835864003

Australia Headquarters

21b Lawler Drive, Oran Park, NSW,
2570, Australia
info@darsunscientific.com
+61-283135514, 430583632

DarsunScientific




Horizontal Electrophoresis

Darsun Scientific

Horizontal Electrophoresis Tank

DS-SHET

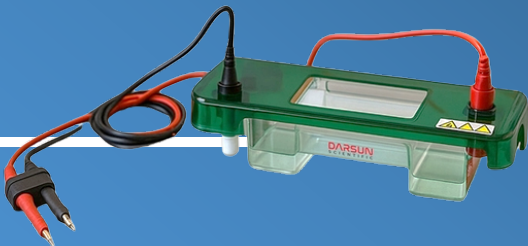

Product Features

Abundant buffer not only guarantees cooling effects, but also keeps the pH value stable during the whole experiment process. Removable electrodes enable maintenance and replacement easy and quick. With background color design on the tray, convenient to adding samples. It's more convenient to open the lid with the easy-open buttons. It prevents buffer volatilizing and electricity leakage with transparent lid. Auto-switch-off when the lid is moved

Model	DS-SHET
Gel size (W x L) No.	48x 75 (mm)
Of samples	3~18
Total buffer volume	~150 (ml)
Dimension (LxWxH)	250x90x50 (mm)
Weight	0.4 (kg)

DS-SHET Packing List

Parts name	Quantity
Lid with cables	1
Body tank with electrodes	1
48x 75mm Gel tray	1
Gel box	1
1.0mm thickness 3-well comb	2
1.5mm thickness 3-well comb	2
1.5mm thickness 5-well comb	2
1.0mm thickness 9-well comb	2
1.5mm thickness 9-well comb	2