




CO₂

INCUBATOR

Darsun Scientific

CO₂ INCUBATOR DS-Co2

DS-50Co2/DS-180Co2/DS-260Co2

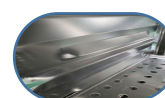

MODEL	DS-50Co2		DS-180Co2	DS-260Co2
Temperature Control Method	Direct heat & Air Jacket			
Temperature Range	RT+3-60°C			
Temperature Uniformity	±0.25°C			
Temperature Accuracy	±0.1°C			
Temperature Recovery Time	<4 Mins			
Ambient Temperature Range	18-34°C			
CO2				
CO2 Range	0 ~ 20%			
CO2 Accuracy	± 0.1%			
CO2 Sensor	Infra-red (IR)			
CO2 Recovery Time*	<4 Mins	<4 Mins	<5 Mins	
Inlet Pressure	0.8-1.2 bar			
HUMIDITY				
Humidification Method	Humidity Pan			
Humidity Range*	>90% RH			
CONTAMINATION				
Contamination method	Dry Heat			
Contamination Temperature	150°C			
Contamination Time	<10h			
Inlet Filter	<0.2µm			
Inner Filter	High temperature resistant HEPA filter, no need to remove during sterilization			
Heater	6 side heating			
Main Body	Antibacterial powder, antibacterial properties >99%			
PHYSICAL CONSTRUCTION				
Interior Volume	50L	185L	260L	
External dims (WxDxH)	500*635*690	685*820*870	745*820*1055	
Internal dims (WxDxH)	310*380*470	500*570*650	555*565*835	
Interior Material	Stainless Steel type 304			
Number of shelves	2	3	4	
Shelves Area	305*270 mm	496*440 mm	550*460 mm	
Maximum Load per Shelf	20kg			
Voltage	220-240VAC 50/60Hz			
Maximum power consumption	500W	1200W	1500W	
Display	Touch screen display with PID controller			
Optional function	i. Co2 Control 1-21%			
	ii. Three inner doors or six inner doors			
	iii. Lot module			
	iv. Humidity Display			
	v. Cylinder Switching			
	vi. Removable base			
	vii. Relief valve			

INTRODUCTION

CONTROL SYSTEM

- Seven inch touch screen, easy to operate, intuitive display
- Three levels of management
- Running record storage, can be exported by USB
- Sound and light alarm, guard sample safety

HEPA FILTER

- High temperature resistant HEPA Filter, no need to remove during sterilization

DOOR SWITCH

- Intelligent door control, adjust the device operation according to the door status

STACKING FEET

- Save lab space

Co2 SENSOR

- I.R sensor dual wave length

SAMPLE PORT

- Easy to sample gas inside the box

DIRECT HEAT & AIR JACKET

- Six sides heating, uniform temperature
- Quick recovery of opening and closing parameters, no overshoot
- 150 or 180° dry contamination

HUMIDITY PAN

- Independent control
- Ensure rapid humidity Recovery

INTERIOR

- Round corner easy to clean, no frame structure

MAIN BODY

- Anti-microbial quoting

SHELVING

- Adjustable steel shelve loading capacity up to 20 kg

STANDING LEG

- Adjustable leg

GLASS DOOR

- For observing sample cells inside the chamber during operation
- protect the clean environment for sample growth





India Office

16 S/F B/S, T-Block Extn.
Jain Colony Part-2 Uttam
Nagar, New Delhi-110059
info@darsunscientific.in
(+91) 9999136670,
7835864003

Australia Headquarters

21b Lawler Drive, Oran Park,
NSW, 2570, Australia
info@darsunscientific.com
+61-283135514, 430583632



DARSUN
S C I E N T I F I C
PRECISE | RELIABLE | ACCURATE

www.darsunscientific.in