



AC

AUTOMATIC CHEMISTRY ANALYZER



India Office

16 S/F B/S, T-Block Extn.
Jain Colony Part-2 Uttam Nagar,
New Delhi-110059
info@darsunscientific.in
(+91) 9999136670, 7835864003

USA Headquarters

30 N Gould St, Ste R
Sheridan, Wyoming 82801
United States

DarsunScientific

System functions

- Throughput: 200 tests /hour
- Automatic, discrete, Randon
- Open the various reagents
- Methods: end-point, two point, kinetic



Reaction system

- 60 pcs reusable plastic cuvettes
- Optical length: 180–400 μl
- Automatic blank testing for cuvettes
- Colorimetry system in the pattern of light splitting



Sample / Reagent / Mixing Handling

- 18 pcs sample position
- 26 pcs reagent posting
- Reagent refrigerated function
- Independent mixing probe
- Sample volume: 1.6–50 μl, 0.1 μl step
- Reagent volume: 10–500 μl, 0.5 μl step
- Auto/manual pre-dilution



Cleaning system

- Auto 7-step washing function
- Low water consumption: 1–2 L/H

Intellective software


- Auto washing cuvette when startup & shut down
- Auto detect reagent volume during every startup
- Reagent plate temperature detecting when startup
- Real-time online monitoring running status reagent
- Residual volume and waste volume

Species

- Dog, cat, horse, monkey, mouse, rabbit
- Pig, cow, goat and sheep

Optical system

- Halogen-tungsten lamp
- 9 wavelengths, 340–691 nm precision ±1.5 nm
- Absorbance range: 0–5.0 abs
- Resolution: 0.0001 Abs rear spectrophotometry

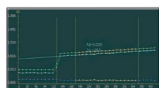

Darsun Scientific

Automatic Chemistry Analyzer

DS-VET

- 200 tests/hour
- Compact and special for VET
- Independent mixing probe
- Minimal labor required
- Smart and easy-operating software



Calibration and auto dilute and re-test

- Linear /nonlinear multi-points calibration
- Auto dilute and retest: – retest when result out of the linearity
- Range or sample is not sufficient, when linear limit is out of range
- Freely set the auto dilute ratio and dilute water position

Input and output

- Input: RS232 interface /computer
- Output: multi-format printout

Operating environment

- Temperature: 10–35 °C
- Humidity: ≤90% no dew

Power requirement

- AC 100–240V
- 50/60 HZ 300VA

Dimension

- 548 mm (H) 422 mm (W) 383 mm (D)

Net weight

- 23 kg