

## Horizontal Electrophoresis Tank

---



### India Office

16 S/F B/S, T-Block Extn.  
Jain Colony Part-2 Uttam Nagar,  
New Delhi-110059  
[info@darsunscientific.in](mailto:info@darsunscientific.in)  
(+91) 9999136670, 7835864003

### USA Headquarters

30 N Gould St, Ste R  
Sheridan, Wyoming 82801  
United States


**Darsun**Scientific

# Horizontal Electrophoresis

Darsun Scientific

## Horizontal Electrophoresis Tank

DS-SHET-A

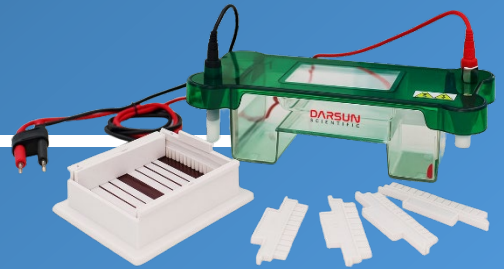

### Product Features

Abundant buffer not only guarantees cooling effects, but also keeps the pH value stable during the whole experiment process.

Removable electrodes enable maintained and replacement easy and quick.

With background color design on the tray, convenient to adding samples.

It's more convenient to open the lid with the easy-open buttons.

It prevents buffer volatilizing and electricity leakage with transparent lid. Auto-switch-off when the lid is moved.

MODEL	DS-SHET-A
Gel size (W x l) No.	48x 75 (mm)
of samples	3-18
Total buffer volume	~150 (ml)
Dimension (LxWxH)	250x90x50 (mm)
Weight	0.4 (kg)

### DS-SHET-A Packing list

Parts name	Quantity
Lid with cables	1
Body tank with electrodes	1
48x 75mm Gel tray	1
Gel box	1
1.0mm thickness 3-well comb	2
1.5mm thickness 3-well comb	2
1.5mm thickness 5-well comb	2
1.0mm thickness 9-well comb	2
1.5mm thickness 9-well comb	2