

Tank Transfer



India Office

16 S/F B/S, T-Block Extn.
Jain Colony Part-2 Uttam Nagar,
New Delhi-110059
info@darsunscientific.in
(+91) 9999136670, 7835864003

USA Headquarters

30 N Gould St, Ste R
Sheridan, Wyoming 82801
United States


DarsunScientific

Horizontal Electrophoresis

Darsun Scientific

Tank Transfer

DS-TB7

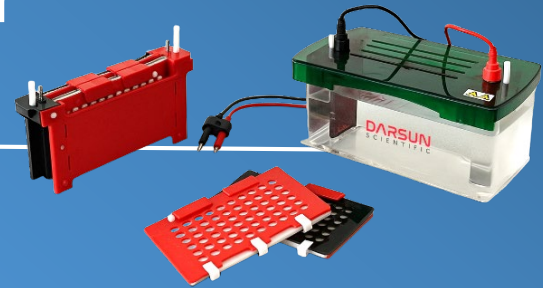

Product Features

Short distance between two electrodes, greatly save transferring time. With color mark on the splint, can ensure to put the gel right into the inner core. Removable electrodes enable maintain and replacement easy and quick. The transparent lid effectively prevent liquid volatilizing and electricity leaking. Auto-switch-off when the lid is moved.

MODEL	DS-TB-7
Trans size (WxL)	216x110(mm)
Total buffer volume	1800 (ml)
Dimension (LxWxH)	297x170x105 (mm)
Weight	1.7 (kg)

Parts name	Quantity
Lid with cables	1
Body tank	1
Core with electrodes	1
Transfer splint	2
Fibre pads	4
Ice box	2