



Touch Screen double-line Printing Sealing Machine with automatic cutter



India Office

16 S/F B/S, T-Block Extn.
Jain Colony Part-2 Uttam Nagar,
New Delhi-110059
info@darsunscientific.in
(+91) 9999136670, 7835864003

USA Headquarters

30 N Gould St, Ste R
Sheridan, Wyoming 82801
United States


DarsunScientific



Darsun Scientific

Touch Screen double-line Printing Sealing Machine with automatic cutter

DS100-CCT

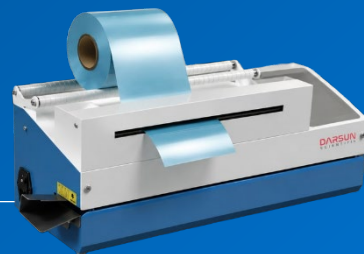

Model	DS100-CCT
Control System	7" color LCD touch screen
Sealing Margin	0-35mm adjustable
Sealing Speed	10±0.5 m/min
Temperature Deviation	≤1%
Cutting Deviation	≤1%
Cutting Width	≤350mm
Sealing Line Width	12mm
Work Temperature	60-220 °C set arbitrarily
Function Type	Automatic cutting, continuous sealing, automatic printing (Double-line print), sealing counting
Heat and Cool	Room temperature rises to 180°C in less than 40s, auxiliary cooling design
Heating Method	Advanced flat heating element, can dry burn, high temperature resistance, long service life, high thermal efficiency
Print Content	Sterilization date, expiration date, pot number, pot batch, operator, items, customized content, etc. Can print English, numbers and international symbols. Sterilization date and expiration date are automatically adjusted and updated every day. 3 adjustable font sizes, printing function can be turned off or a certain item can be selectively closed.
Safety Guarantee	The deviation between the sealing temperature and the set working temperature exceeds 4 °C, it will automatically stop working
Indentation Style	Flat wheel sealing (Default) / Single indentation sealing / Pattern wheel sealing
Pressing System	The floating constant pressure pressing system is designed to meet the sealing needs of various bags